

CEDARON CONNECT


Connect – Steps to Configure Faxing Patient Reports

In order to use the functionality available in Connect for faxing patient reports, the following are required:


1. A fax modem
2. Installing Windows Fax Services on the client pc (that is done in windows add/remove windows components)
3. A "Fax Mail Transport" account to be set up in Microsoft Outlook

Setting up a Fax Mail Transport account:

Open Outlook and add an Email Account

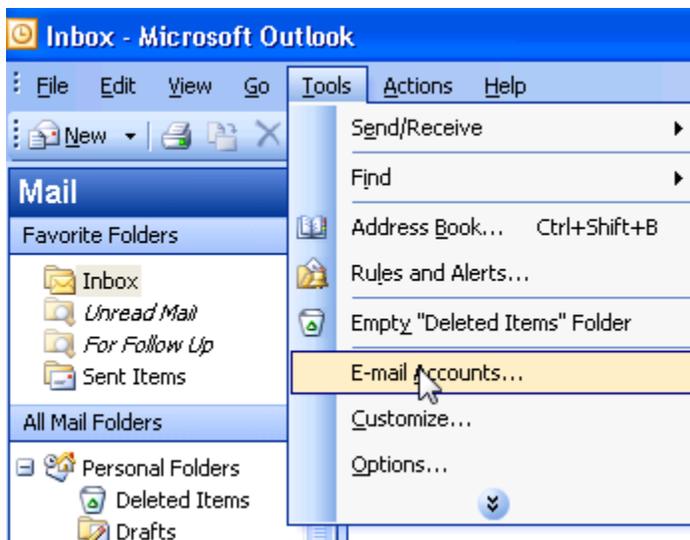

CEDARON CONNECT

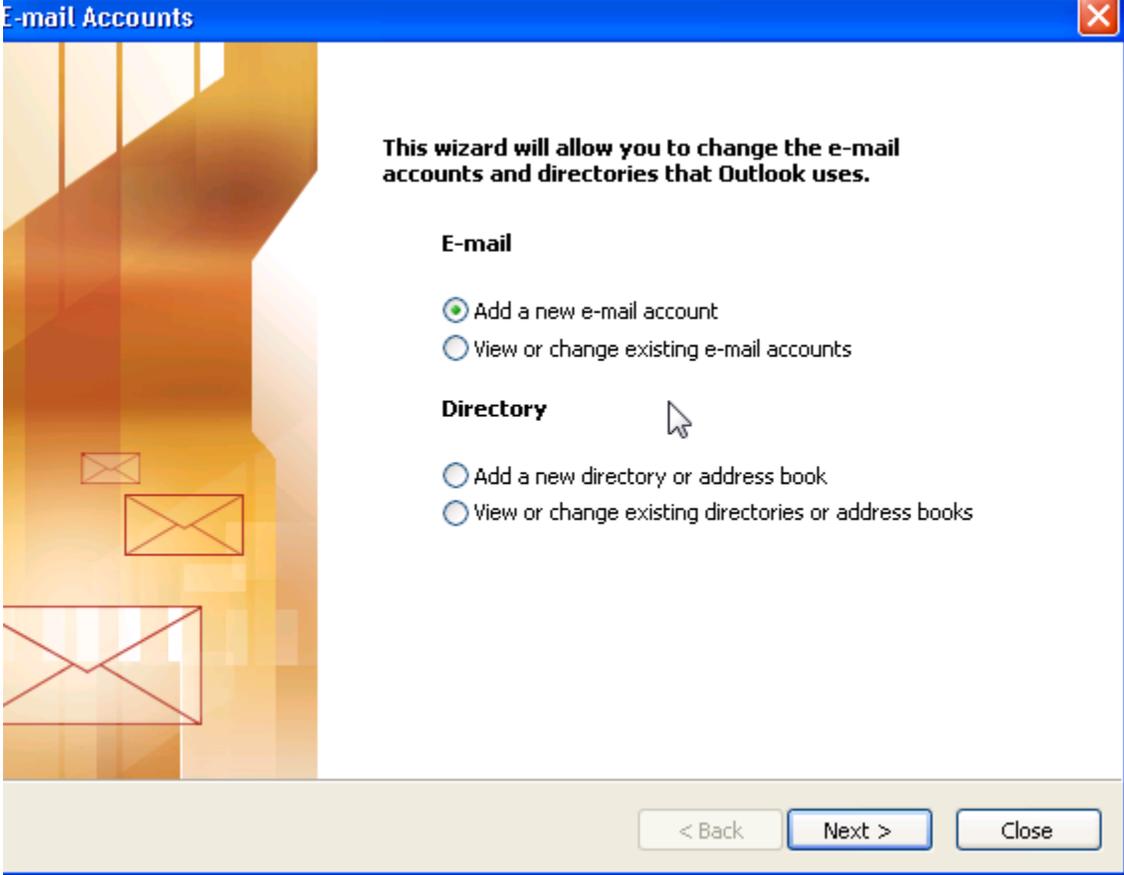

CEDARON CONNECT

Choose "Additional Server Types"

E-mail Accounts [Close]


Server Type
You can choose the type of server your new e-mail account will work with.

- Microsoft Exchange Server**
Connect to an Exchange server to read e-mail, access public folders, and share documents.
- POP3**
Connect to a POP3 e-mail server to download your e-mail.
- IMAP**
Connect to an IMAP e-mail server to download e-mail and synchronize mailbox folders.
- HTTP**
Connect to an HTTP e-mail server such as Hotmail to download e-mail and synchronize mailbox folders.
- Additional Server Types**
Connect to another workgroup or 3rd-party mail server.

< Back **Next** Cancel

CEDARON CONNECT

Select "Fax Mail Transport"

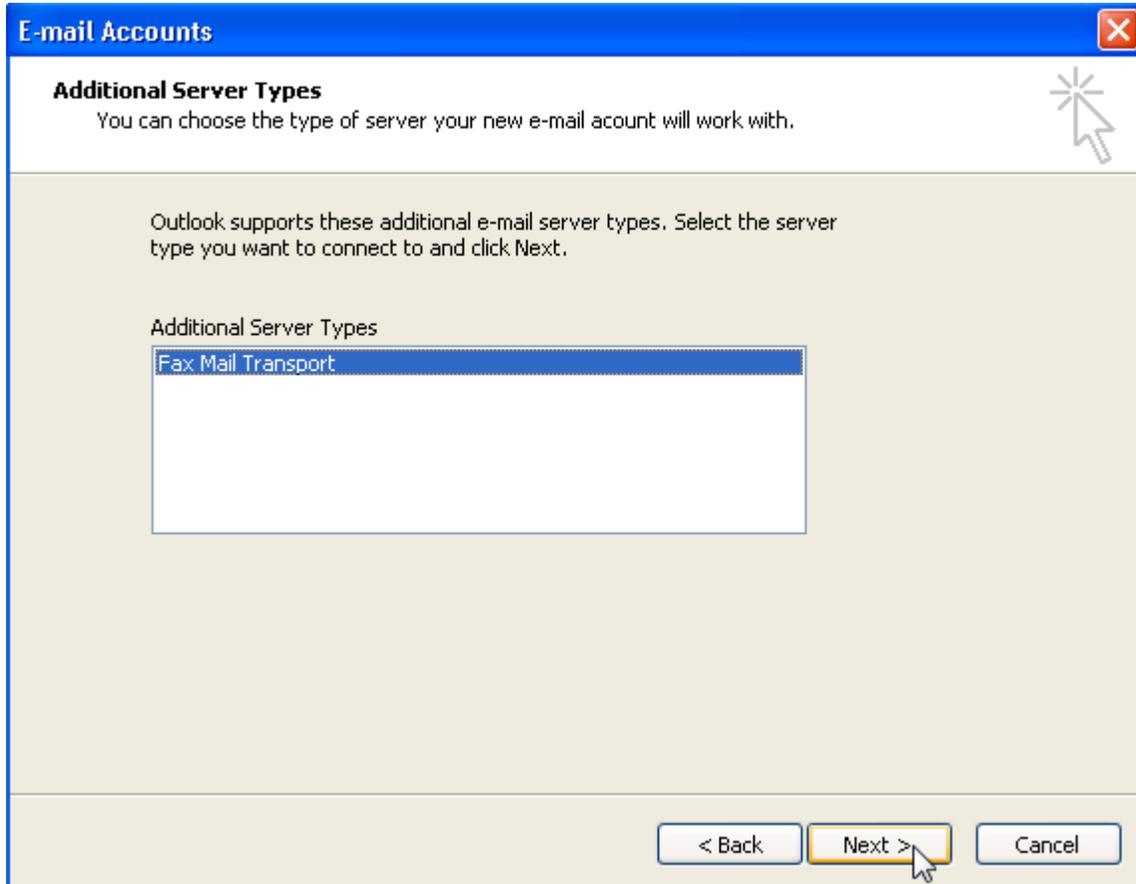

Close Outlook.

If you are using Microsoft XP you will need to download and install Windows Search 4.0


Create the following registry key ([HKEY_LOCAL_MACHINE\SOFTWARE\Microsoft\At Work Fax]) by going to Start -> Run -> regedit

CEDARON CONNECT

Expand HKEY_LOCAL_MACHINE

Expand SOFTWARE

Right click on Microsoft on select New -> Key

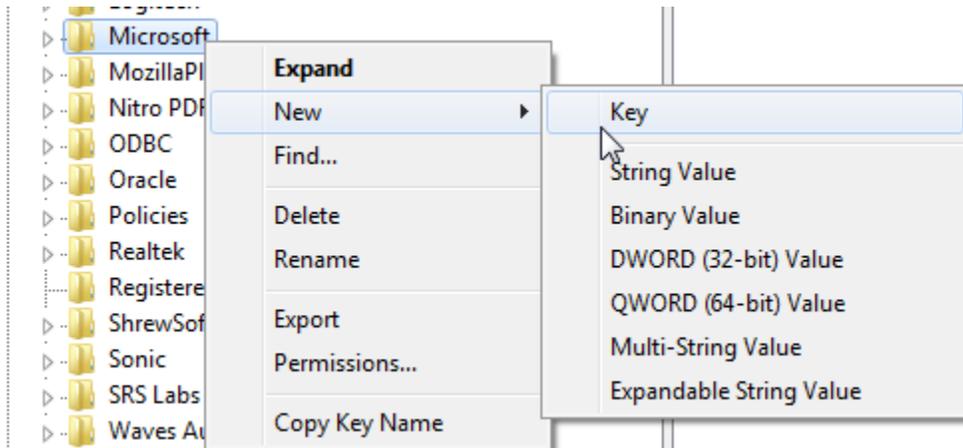

The new key name should be "At Work Fax" (no quotes).



Now you can start Outlook and send a test fax if desired. Connect should be ready to fax now.